



Silk Line™

SOFT DENTURE RELINING MATERIAL

Silk Line

Silk Line is a long-lasting, soft denture relining material based on addition-curing vinyl silicones. Silk Line is indicated for use with all polymethylmethacrylate-based synthetics. Delivered in an auto-mix cartridge, Silk Line allows quick and easy chairside soft relines, in addition to being indispensable in the dental laboratory.

Silk Line relieves pressure areas and sore spots in the upper and lower jaws. It also improves retention by creating a functional post-dam and adapts well even to flabby ridges. At the same time, sharp-ridged alveolar processes are elastically covered, and increased retention and stability improve the occlusion of full upper and lower dentures. Used in conjunction with overdentures during the healing process, Silk Line is also suitable for cushioning implants.



Features and Benefits

- Auto-mix cartridge permits direct application to denture base
- Unique primer enhances durability
- Composition allows easy trimming and polishing
- Relines can withstand wear up to two years
- Chairside use is convenient for patients
- Absence of odor and taste increases patient satisfaction
- May also be used in a dental laboratory

**To order, or for more information,
call 888-JMORITA (566-7482) or contact
your dealer.**

J.MORITA
USA INC.

Regional Partner


MORITA

Silk Line Primer Improves Durability

The difficulty of adhering a silicone-based product to an acrylic denture base is eliminated with Silk Line's primer.

This unique material was formulated to create a resilient bond between the PMMA chemistry of the denture and the silicone contained in the Silk Line paste. The result is a very durable, yet soft reline, resistant to long term use and cleaning.

Varnish (Optional Use)


The Silk Line Kit contains a varnish catalyst and a varnish base. The varnish can be used in the final polishing step to smooth the trimmed borders of the denture, if desired. The base and catalyst are mixed together in a 1:1 ratio, brushed onto the denture, and allowed to dry at room temperature for approximately 5 minutes.



Dispensing gun sold separately

Basic Steps For Use of Silk Line:

- 1) Check fit of denture
- 2) Remove 1-2 mm from its acrylic base
- 3) Apply primer
- 4) Dispense Silk Line paste onto denture
- 5) Insert denture into patient's mouth
- 6) Remove after 4 minutes
- 7) Finish and polish

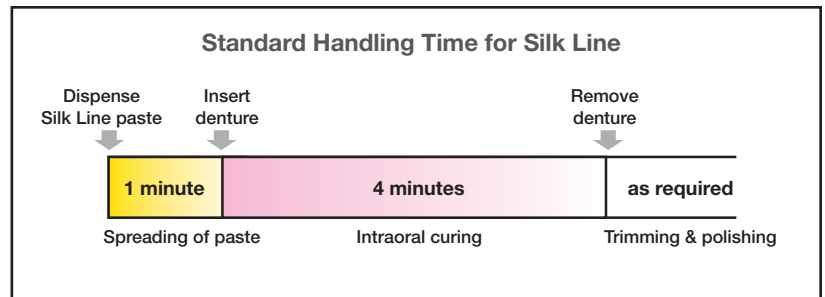

Ordering Information

Packaging

Kit contains: Auto-Mix Cartridge (50ml), Primer (5ml), Varnish Catalyst (5ml), Varnish Base (5ml), Mixing Tips (10), Disposable Mixing Well (10), Disposable Brushes (20), and Brush Handle (1)

Product Numbers

- Silk Line Kit.....28-100100
- Silk Line Refill Cartridge (50ml)28-100200
- Silk Line Primer (5ml)28-100300
- Mixing Tips Refill (12)21-300047
- Mixing Tips Refill (50)21-300047L

