

SIGNO  
GIOII



Thinking ahead. Focused on life.



# Signo G10 II

## **The evolution has begun**

The new Signo G10 II is the latest evolution of the proven Signo range of dental treatment units from Morita. Signo G10 II offers a contemporary streamlined aesthetic appeal, complemented by a range of advanced dental elements designed to meet the requirements of today's dental professional.



### **Luna Vue ES Operating Light**

Signo G10 II offers a new integrated LED operating light which provides vivid natural colour with high luminescence, while still providing practical advantages of high energy efficiency and long lasting durability.



### **Chair Design**

In addition to contemporary aesthetic appeal, the Signo G10 II chair has been designed to provide maximum flexibility for patient access. The height of the chair can be adjusted from a minimum of 400mm to a maximum of 750mm. This wide range ensures easy access for patients of all ages.



### **Tray Arm**

The operator's tray and instrument holder support arm is designed to provide maximum flexibility to allow easy patient access to the treatment area. The operator can easily move the tray vertically and horizontally to match their desired operating position. While an arm movement lock ensures a safe environment for both patient and operator.

### New Continental Delivery Option

Signo G10 II's continental delivery highlights 5 instrument lines with perfectly balanced, self retracting continental style handpiece levers. The newly designed straight-handled triplex syringe seamlessly adapts to the whip arm delivery. Optional accessory trays are also available to complement the new continental delivery system to optimize workflow performance.



### Instrument Holder and Tray

The Signo G10 II instrument holder provides a seamless transition to the generously sized operator's tray. The instrument holder incorporates a flexible array of instruments including a three-way syringe, highspeed handpieces, an electric micromotor and scaler. Optionally available with TwinPower or TorqTech handpieces.



# Basin

## Glass Basin

Focused on quality and cleanliness the Signo G10 II offers a contemporary glass basin design. The basin is easily removed for cleaning and is available in a range of colours to complement your upholstery choice.



Clear



Green



Orange



Pink

# Upholstery and Trim

## Customized Appeal

The Signo G10 II incorporates further scope to customize the appearance of your treatment unit with 2 new trim options (tray and base sections). This combined with multiple upholstery and basin colour options ensures your Signo G10 II will fit perfectly into your practice with maximum appeal.



G10 Gloss Black (G10 GB)



G10 Shiny Silver (G10 SS)



G10 Fresh Green (G10 FG)



G10 Cool Violet (G10 CV)



G10 Sweet Brown (G10 SB)



G10 Splash Yellow (G10 SY)



G10 Honey Pink (G10 HP)

# Standard

# Options



Three-way syringe (x2)  
(WS-201)



Highspeed handpiece  
TwinPower (with optic)  
TwinPower Ultra Series (with optic)



Luna Vue EL II Light  
(12 LED)



Small tray mounted on light pole  
Tissue box holder  
(Over arm type)



Vacuum syringe (VS14WD2)  
(with suction power control valve)



Electric micromotor  
TORX TR-91-O (optic)



EL mirror



Foot controller  
(Extended function with lever)  
(FC201A)



Luna Vue ES Light  
(Single LED)



AIR TORX TRA-200



Movable armrest  
Dentist's and Assistant's side




Triple viewer (PFV95LD)

### Additional standard inclusions:

- 2 x Highspeed (Midwest) handpiece lines
- Electric micromotor (TR-91)
- Standard foot controller
- Saliva ejector

Specifications	
Name	SIGNO G10 II
Model	TU101
Input voltage	Adjusted to local line voltage
Power consumption	1.5kVA
Patient seat elevation system	Hydraulic motor
Seat height	750mm max. 400mm min.
Backrest angle (from horizontal line)	Raised 70° Lowered 2°
Weight	Approx. 240kg

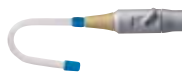

Saliva ejector  
SE14S  
(with suction power control valve)



Armrest  
Assistant's side



Operating stool LN II (OS-LN II)  
Width: 510mm (seat width 400mm)  
Depth: 540mm (seat width 450mm)  
Seat height: 445-635mm  
Height adjustable range: 190mm  
Weight: Approx. 16kg



Articulated headrest



Bacteria filters  
(water and air)

### Additional options:

- Monitor mount (light pole)
- Accessory trays
- Water line flushing system

Diagnostic and Imaging Equipment

Treatment Units

Handpieces and Instruments

Endodontic Systems

Laser Equipment

Laboratory Devices

Educational and Training Systems

Auxiliaries



Development and Manufacturing

**J. Morita Tokyo Mfg. Corporation**

7129 Komuro, Ina-machi, Kitaadachi-gun,  
Saitama 362-0806, Japan

T +81. (0)48.723 2621, F +81. (0)48.723 2622

**Morita Global Website**

**[www.morita.com](http://www.morita.com)**

Distribution

**J. Morita Corporation**

3-33-18 Tarumi-cho, Suita-shi, Osaka 564-8650, Japan

T +81. (0)6. 6380 1521, F +81. (0)6. 6380 0585

**J. Morita Europe GmbH**

Justus-von-Liebig-Strasse 27a, 63128 Dietzenbach, Germany

T +49. (0)6074. 836 0, F +49. (0)6074. 836 299

**Morita Dental Asia Pte.Ltd.**

3 Science Park Drive,

#01-05 The Franklin Singapore Science Park1,

Singapore 118223

T+65. 6779. 4795, F+65. 6777. 2279

**J. Morita Corporation Australia & New Zealand**

Suite 2.05, 247 Coward Street, Mascot NSW 2020, Australia

T +61. (0)2. 9667 3555, F +61. (0)2. 9667 3577

**J. Morita Corporation Middle East**

4 Tag Al Roasaa, Apartment 902, Saba Pacha 21311 Alexandria, Egypt

T +20. (0)3. 58 222 94, F +20. (0)3. 58 222 96

**J. Morita Corporation India**

Felix Office No.908, L.B.S. Marg, Opp. Asian Paints, Bhandup (West), Mumbai 400078, India

T +91. 22. 2595. 3482

**J. Morita Mfg. Corporation Indonesia Representative Office**

28F, Ciputra World Tower 1, Jl. Prof. Dr. Satrio Kav. 3-5, Jakarta 12940, Indonesia

T +62. 21. 2988 8287, F +62. 21. 2988 8201

**Siamdent Co., Ltd.**

444 Olympia Thai Tower, 3rd Floor, Ratchadapisek Road.

Samsennok, Huay Kwang, Bangkok 10310, Thailand

T +66. (0)2. 512 6049, F +66. (0)2. 512 6099

[www.siamdent.com](http://www.siamdent.com)